

MiRo House Becklands Close Boroughbridge North Yorkshire YO51 9NR UNITED KINGDOM

tel: +44(0)1423 320 007

# TOR MAS / MAX Mammography Phantoms

A test object designed to be used quickly and easily on a routine basis to provide an ongoing check of imaging performance, particularly those aspects which are most liable to deterioration. An ongoing record of these numbers will reveal any trend towards deterioration in imaging performance.

### TOR MAS and TOR MAX Consist of the following:

#### Standard PMMA test block

(6 pieces 10.0±10%mm thick pmma 240mm diameter semi-circle 2 pieces 5.0±10%mm thick pmma 240mm diameter semi-circle)

### • Sensitometry

(Ten-step grey-scale plus two points for Sensitometric measurements)

### High Resolution limit

(1.0 to 20.0 LP/mm) x1 for TOR MAS x2 for TOR MAX

#### Low Contrast Resolution

(1.8 to 5 line pairs/mm, representing filamentary structures)

## Low-contrast large-detail detectability

(12 details, 5.6mm diameter)

### • High-contrast small-detail detectability

(11 details, 0.5 and 0.25mm diameter)

#### Microcalcifications

(Particles representing microcalcifications on a step-wedge background, median sizes 125, 234, 328 microns)



深圳市展业达鸿科技有限公司

吕先生:15920060912(微信同号) 0755-22934005(座机)

地址:深圳市福田区八卦二路八卦岭工业区615栋419

邮箱:hongqi@thingstet.com




MiRo House Becklands Close Boroughbridge North Yorkshire YO51 9NR UNITED KINGDOM


tel: +44(0)1423 320 007

## TOR MAS / MAX Mammography

# **Product X-Rays**

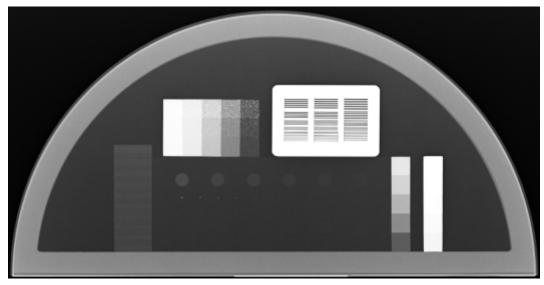

### fig. 1 TOR MAS X-ray image

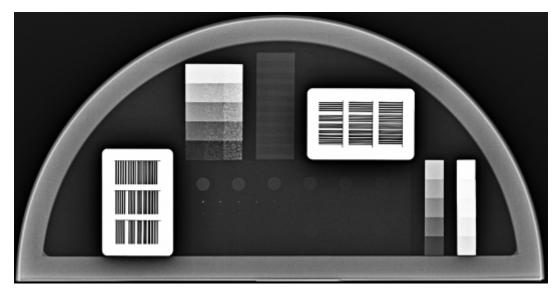

### fig. 2 TOR MAX X-ray image

深圳市展业达鸿科技有限公司 日先生:15920060912(微信同号) 0755-22934005(座机)

地址:深圳市福田区八卦二路八卦岭工业区615栋419

邮箱:hongqi@thingstet.com

Copyright © 2019 Leeds Test Objects Ltd | SpecSh 179 题物的/2009 thingstest.com